

Ureteral Catheters



Used for drainage of urine from the ureter and it also allows the delivery of contrast media.

- Ureteral Catheters are radiopaque with centimeter marking
- Ureteral Catheters are supplied with ureteral catheter adapter
- Angled Open Tip & Angled Closed Tip supplied with Stylet

Open Tip

- Open tip allows placement over a guidewire for negotiating a tortuous ureter

Open Tapered Tip

- Open tip allows placement over a guidewire for negotiating a tortuous ureter
- Tip is tapered for ease of insertion into ureteral orifice

Closed Tip

- Closed tip allows placement over a guidewire for negotiating a tortuous ureter
- Smooth rounded tip
- Stylet allows negotiating through tortuous ureter

Cone Tip

- The Tip occludes Ureteral orifice for retrograde pyelogram

Angled Open Tip

- Open tip allows placement over a guidewire for negotiating a tortuous ureter
- Angle Tip allows manipulation during insertion

Angled Closed Tip

- Angled Tip allows manipulation during insertion
- Stylet allows negotiation through tortuous ureter



Open Tip




Open Tapered Tip




Closed Tip



Cone Tip

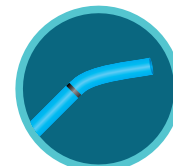

Angled Open Tip

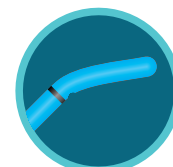

Angled Closed Tip

QTY. / Box : 10

Order No. Open Tip	Order No. Open Tapered tip	Order No. Closed Tip	Order No. Cone tip	Order No. Angled Open Tip	Order No. Angled Closed Tip	Fr	Length (cm)
UC-O3	-	UC-R3	-	UC-AO3	UC-AR3	3	70
UC-O4	-	UC-R4	-	UC-AO4	UC-AR4	4	70
UC-O5	UC-OT-5	UC-R5	UC-C5/8	UC-AO5	UC-AR5	5	70
UC-O6	UC-OT-6	UC-R6	UC-C6/10	UC-AO6	UC-AR6	6	70
UC-O7	UC-OT-7	UC-R7	UC-C7/12	UC-AO7	UC-AR7	7	70